

Specimen Instructions

Peripheral Whole Blood


Use only tubes provided inside the FoundationOne®Liquid Collection and Shipping Box. Other tubes will not be accepted.



Overview of Specimen Handling

Accurate analysis of cell-free DNA requires proper collection technique and handling of the sample. Failure to adhere to these instructions can compromise results by diluting cell-free DNA with DNA from white blood cell lysis.

- 1** Check special tubes provided in FoundationOne®Liquid Collection and Shipping Box to confirm liquid is clear and without cloudiness or crystals.
- 2** Label each tube with the supplied labels indicating the date of collection and two unique patient identifiers (order identification number and date of birth).
- 3** Collect two tubes of whole blood (8.5mL per tube).
 - Prevent backflow: tubes contain chemical additives and it is important to avoid backflow into patient.
 - Collect specimen by venipuncture according to CLSI H3-A6.1
 - Fill tubes to marked level (8.5mL per tube).
- 4** Remove the tube from adapter and immediately **mix by gentle inversion 8 to 10 times**. Inadequate or delayed mixing may result in inaccurate test results. One inversion is a complete turn of the wrist, 180 degrees, and back per the figure below.

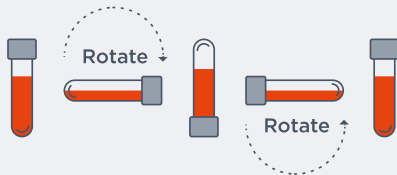

- 5** Place specimen, matching Test Requisition Form (TRF), and Diagnosis Confirmation sheet into the FoundationOne®Liquid Collection and Shipping Box
 - Confirm each tube is labeled with the supplied labels indicating the date of collection and two unique patient identifiers (order identification number and date of birth).
- 6** Preferably on the same day of collection, arrange for courier pick-up. Do not freeze or refrigerate blood samples.

Temperature is important. Keep at room temperature (18-25 degrees Celsius).
DO NOT FREEZE.

Contact courier to send sample to:

Foundation Medicine®, Inc.
150 Second Street
Cambridge, MA 02141
Phone: +49 (0)7624 / 14-2098
Email: europa.foundationmedicine@roche.com

* See separate Packing Instructions for details.

Reference

1. Clinical and Laboratory Standards Institute. H3-A6, Procedures for the collection of diagnostic blood specimens by venipuncture; approved standard-sixth edition.

